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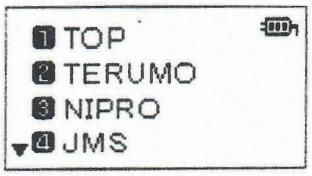
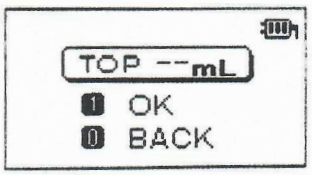
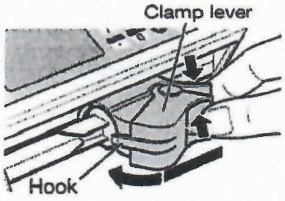
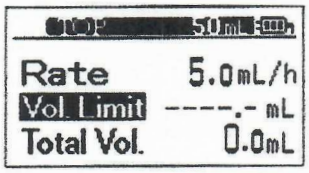
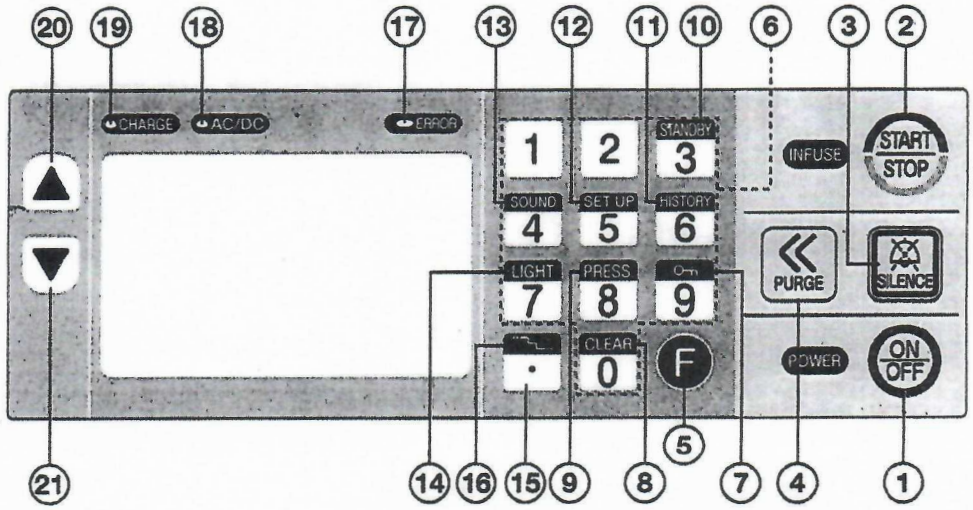
No.	Name	Function	No.	Name	Function
①*1	Power key	Turns power supply on/off.	⑫*2	Setting key	Enters the Setting Mode.
②	Start/Stop key	Starts/Stops the infusion.	⑬*2	Volume key	Displays a screen for adjusting the buzzer volume.
③	Silence key	Silences the alarm.	⑭*2	Light key	Turns the backlight on/off.
④	Purge key	Performs purging.	⑮	Decimal point key	Use to set a decimal point value.
⑤	Function key	Used for special functions.	⑯*2	Infusion rate change key	Used to change the flow rate while infusion is under way.
⑥	Numeric key	Used to enter a numerical value.	⑰	Error indicator	Indicates a system fault.
⑦*2	Key Lock key	Locks the operation key.	⑱	AC/DC power supply indicator	Illuminated when external power is supplied.
⑧*2	Clear key	Clears the numerical value.	⑲	Battery charge indicator	Illuminated while the internal battery is being recharged.
⑨*2	Pressure Display key	Displays the pressure state during infusion on a graph.	⑳	UP key	Switches between items when setting numerical values.
⑩*2	Standby key	Makes transition to Standby mode.	㉑	DOWN key	Switches between items when setting numerical values.
⑪*2	History key	Displays the history.			